
TAYLOR®

Helpful Hints



Model# 5345339

IMPORTANT, RETAIN FOR FUTURE REFERENCE: READ CAREFULLY

DOWNLOAD APP

Scan the QR code to get the App



PRECISION HUB
BY TAYLOR

**Download the
Precision Hub App
Scan Code**



To get the latest and
most up-to-date
version of the App
use the QR code.

NOTE: Due to
improvements in software
the APP may change in its
appearance.

09:36

< User Info

Nickname
Please enter Nickname 0/25

Gender

Birthday 06-07-2000 >

Height 5' 2" ft >

Current Weight 150.6lb >

Goal Weight 125.0lb >

Save

Create your profile

09:39

< Devices

Name: Digital Scale
Model: 5345339
Mac: 08:B8:D0:23:42:09 >

+ Click here to add a device.

**Using connecting
technology pair your
phone to the scale**

HOW TO USE YOUR SCALE

1 INSERT THE BATTERIES

- **2** 'AAA' (1.5V) batteries are included with your scale.



2 YOUR SCALE IS SET AT THE FACTORY TO WEIGH IN POUNDS (LB)

- While the scale is on, press the “UNIT” unit button on the back of the scale to convert the weight unit to kg if needed.



3 SET THE SCALE IN THE RIGHT POSITION

- Use only on a flat and hard floor surface. Do not use on carpets.



Tap scale to activate

4 START-UP THE SCALE

- Press firmly once on the scale platform to start up the scale.
- The display shows “0.0” and then turns off.
- The scale is ready for use.


Note: If you move the scale at all, please repeat the steps above for best accuracy.




Step evenly on the platform with bare feet

5 WEIGHING YOURSELF

- With bare feet, step onto the scale and stand still while your weight is being measured.
- The scale will count up to your weight and lock it in. The scale will hold the display with the weight showing.
- The scale turns off automatically after use.

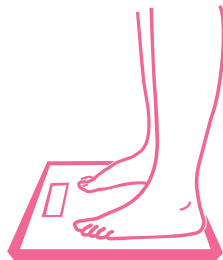

Turn on your connecting technology phone settings then start the app and go to the overview screen. With bare feet, stand on your scale. Continue to stand still on the scale in bare feet.

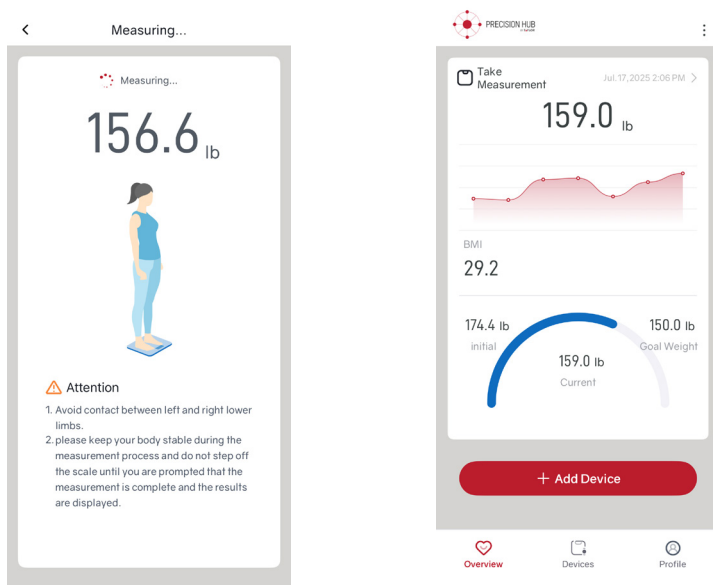

Once your weight and BMI have been measured, they will be displayed in the app. The app will track your changes over time.

Note: The Taylor Smart Scale can be used with and without the Precision Hub App. Also, note that due to constant improvements to the app, the displays may appear different from what is shown in this guide.

TROUBLESHOOTING

If you are having trouble synchronizing data, device connections, or other scale use questions, please refer to the help button in the app or the scale instruction manual for more detailed information.

SAFETY WARNING

- The warning signs and icons shown appear in order for you to use this product safely and correctly, as well as to prevent product damage, risk, and injury to you and others.

PRECAUTION NOTICE

- Do not disassemble the scale, as incorrect handling may cause injury.
- Do not step on the scale when your body or feet are wet, especially after bathing or showering, to prevent slipping.

IMPORTANT NOTICE TO USERS

- Make sure to use only the type of battery stated in the manual.

TAYLOR

2220 Entrada Del Sol, Suite A
Las Cruces, New Mexico 88001 USA

Phone: 1-866-843-3905

Email: taylorusa@lifetimebrands.com

©2026 Taylor, a division of Lifetime Brands, Inc. All rights reserved. Oak Brook, IL,
60523

**Made to our exact specifications in China.
Not legal for trade.**

