

Client Samsung (Dymstec)

CALIBRATION CERTIFICATE

Object(s) ET3DV6 - SN 1617

Calibration procedure(s) QA CAL-01.v2
Calibration procedure for dosimetric E-field probes

Calibration date: September 18, 2003

Condition of the calibrated item In Tolerance (according to the specific calibration document)

This calibration statement documents traceability of M&TE used in the calibration procedures and conformity of the procedures with the ISO/IEC 17025 international standard.

All calibrations have been conducted in the closed laboratory facility: environment temperature 22 +/- 2 degrees Celsius and humidity < 75%.

Calibration Equipment used (M&TE critical for calibration)

Model Type	ID #	Cal Date (Calibrated by, Certificate No.)	Scheduled Calibration
Power meter EPM E4419B	GB41293874	2-Apr-03 (METAS, No 252-0250)	Apr-04
Power sensor E4412A	MY41495277	2-Apr-03 (METAS, No 252-0250)	Apr-04
Reference 20 dB Attenuator	SN: 5086 (20b)	3-Apr-03 (METAS No. 251-0340)	Apr-04
Fluke Process Calibrator Type 702	SN: 6295803	8-Sep-03 (Sintrel SCS No. E-030020)	Sep-04
Power sensor HP 8481A	MY41092180	18-Sep-02 (Agilent, No. 20020918)	In house check: Oct 03
RF generator HP 8684C	US3642U01700	4-Aug-99 (SPEAG, in house check Aug-02)	In house check: Aug-05
Network Analyzer HP 8753E	US37390585	18-Oct-01 (Agilent, No. 24BR1033101)	In house check: Oct 03

	Name	Function	Signature
Calibrated by:	Katja Pokovic	Laboratory Director	
Approved by:	Fin Bortholt	R&D Director	

Date issued: September 18, 2003

This calibration certificate is issued as an intermediate solution until the accreditation process (based on ISO/IEC 17025 International Standard) for Calibration Laboratory of Schmid & Partner Engineering AG is completed.

Probe ET3DV6

SN:1617

Manufactured:	January 25, 2002
Last calibration:	July 16, 2002
Recalibrated:	September 18, 2003

Calibrated for DASY Systems

(Note: non-compatible with DASY2 system!)

DASY - Parameters of Probe: ET3DV6 SN:1617

Sensitivity in Free Space

NormX	1.78 $\mu\text{V}/(\text{V}/\text{m})^2$
NormY	1.75 $\mu\text{V}/(\text{V}/\text{m})^2$
NormZ	1.79 $\mu\text{V}/(\text{V}/\text{m})^2$

Diode Compression

DCP X	96	mV
DCP Y	96	mV
DCP Z	96	mV

Sensitivity in Tissue Simulating Liquid

Head 900 MHz $\epsilon_r = 41.5 \pm 5\%$ $\sigma = 0.97 \pm 5\%$ mho/m

Valid for f=800-1000 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	6.3 $\pm 9.5\%$ (k=2)	Boundary effect:
ConvF Y	6.3 $\pm 9.5\%$ (k=2)	Alpha 0.40
ConvF Z	6.3 $\pm 9.5\%$ (k=2)	Depth 2.60

Head 1800 MHz $\epsilon_r = 40.0 \pm 5\%$ $\sigma = 1.40 \pm 5\%$ mho/m

Valid for f=1710-1910 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	5.3 $\pm 9.5\%$ (k=2)	Boundary effect:
ConvF Y	5.3 $\pm 9.5\%$ (k=2)	Alpha 0.52
ConvF Z	5.3 $\pm 9.5\%$ (k=2)	Depth 2.54

Boundary Effect

Head 900 MHz Typical SAR gradient: 5 % per mm





Probe Tip to Boundary		1 mm	2 mm
SAR _{be} [%]	Without Correction Algorithm	11.1	6.4
SAR _{be} [%]	With Correction Algorithm	0.3	0.6



Head 1800 MHz Typical SAR gradient: 10 % per mm

Probe Tip to Boundary		1 mm	2 mm
SAR _{be} [%]	Without Correction Algorithm	13.5	9.0
SAR _{be} [%]	With Correction Algorithm	0.2	0.1


Sensor Offset

Probe Tip to Sensor Center	2.7	mm
Optical Surface Detection	1.4 $\pm$ 0.2	mm

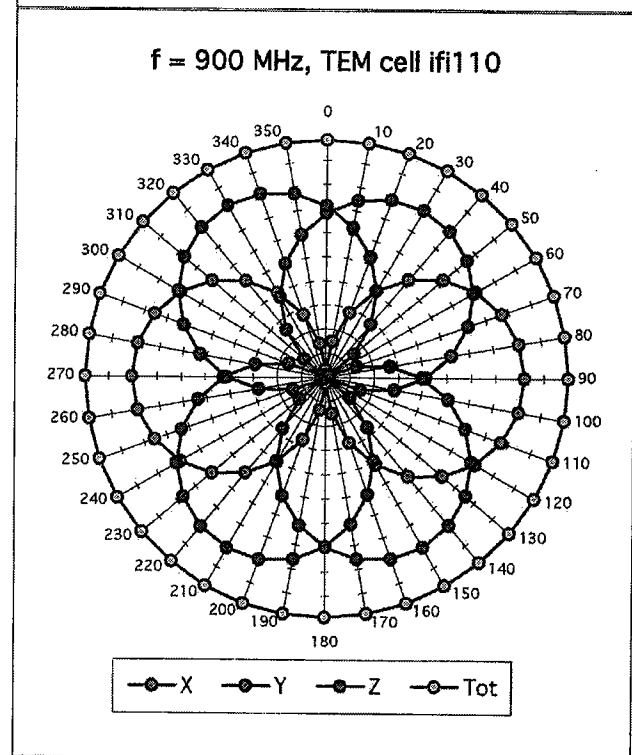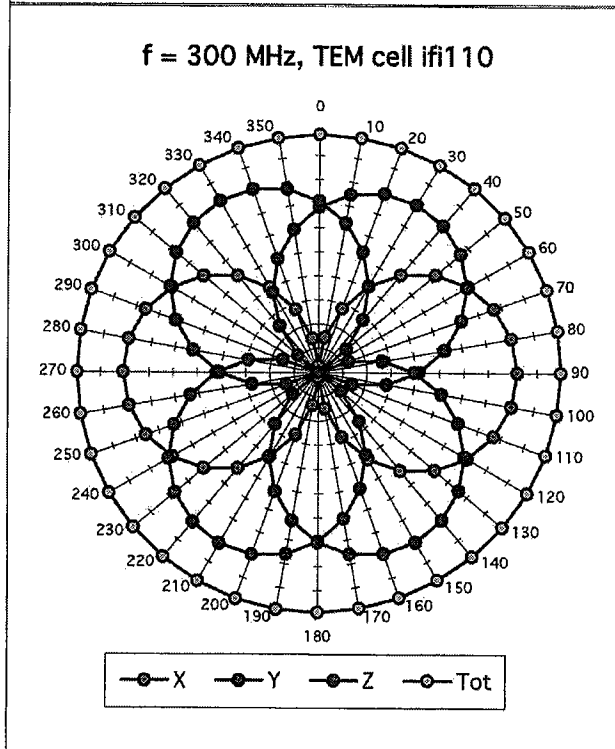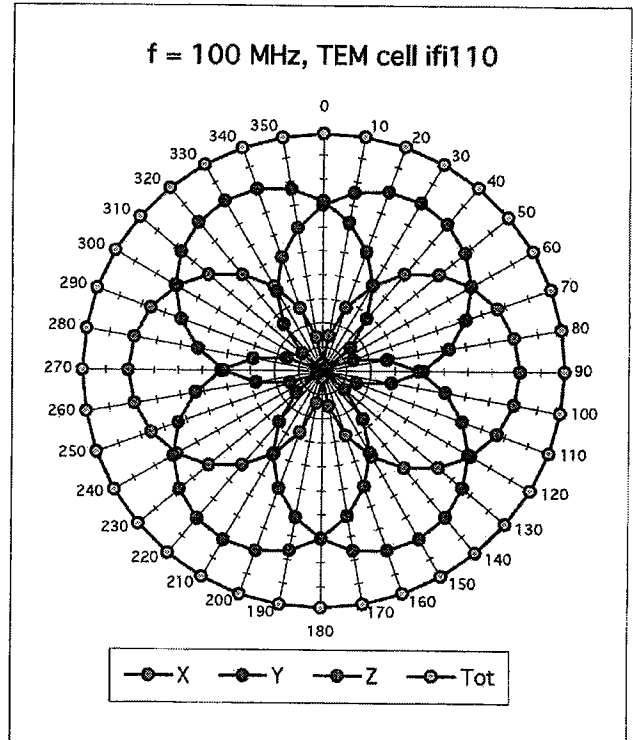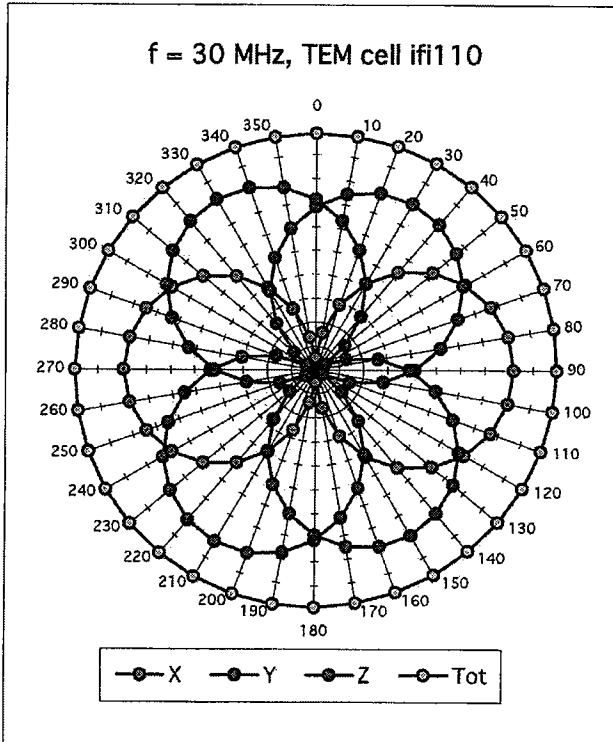Receiving Pattern (ϕ), $\theta = 0^\circ$ 

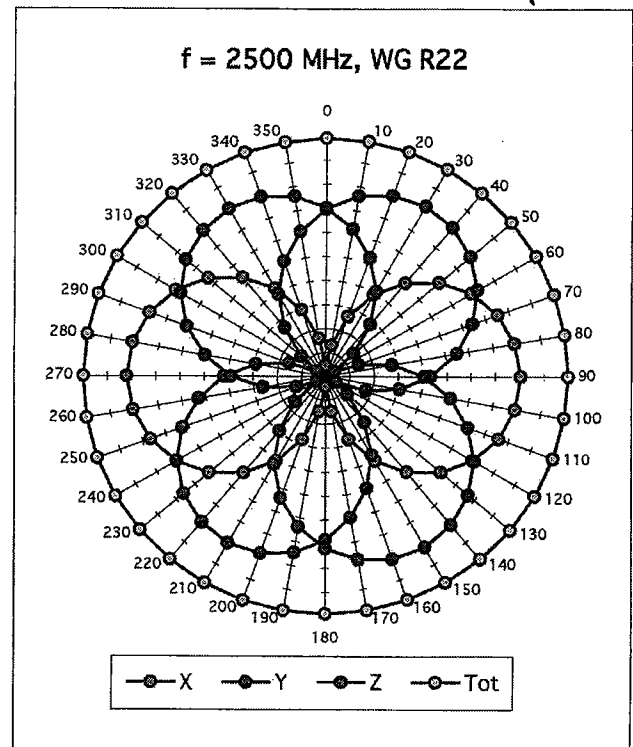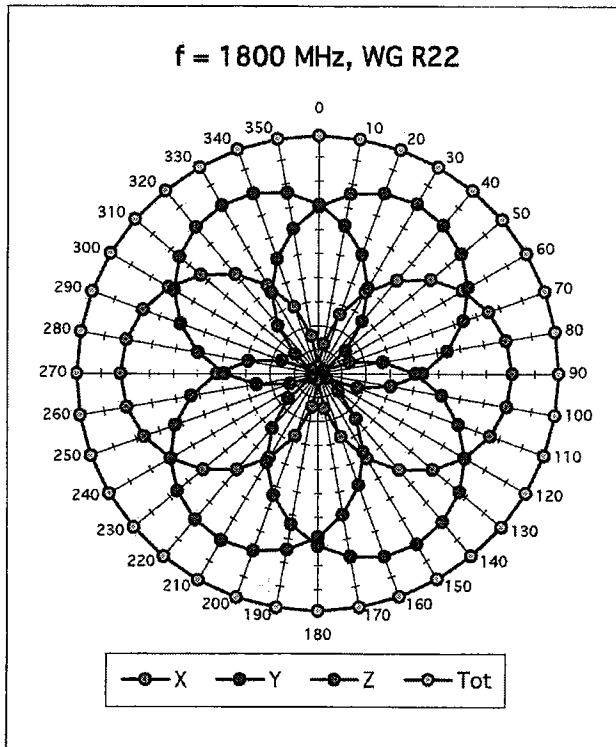


Isotropy Error (ϕ), $\theta = 0^\circ$

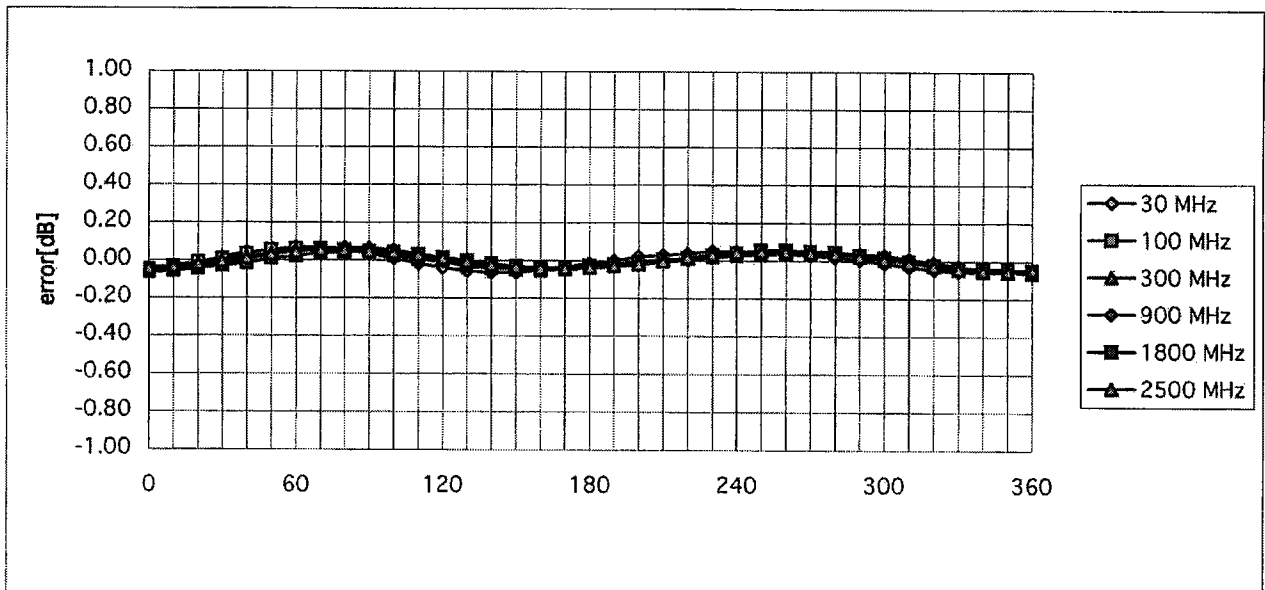



Frequency Response of E-Field

(TEM-Cell:ifi110, Waveguide R22)


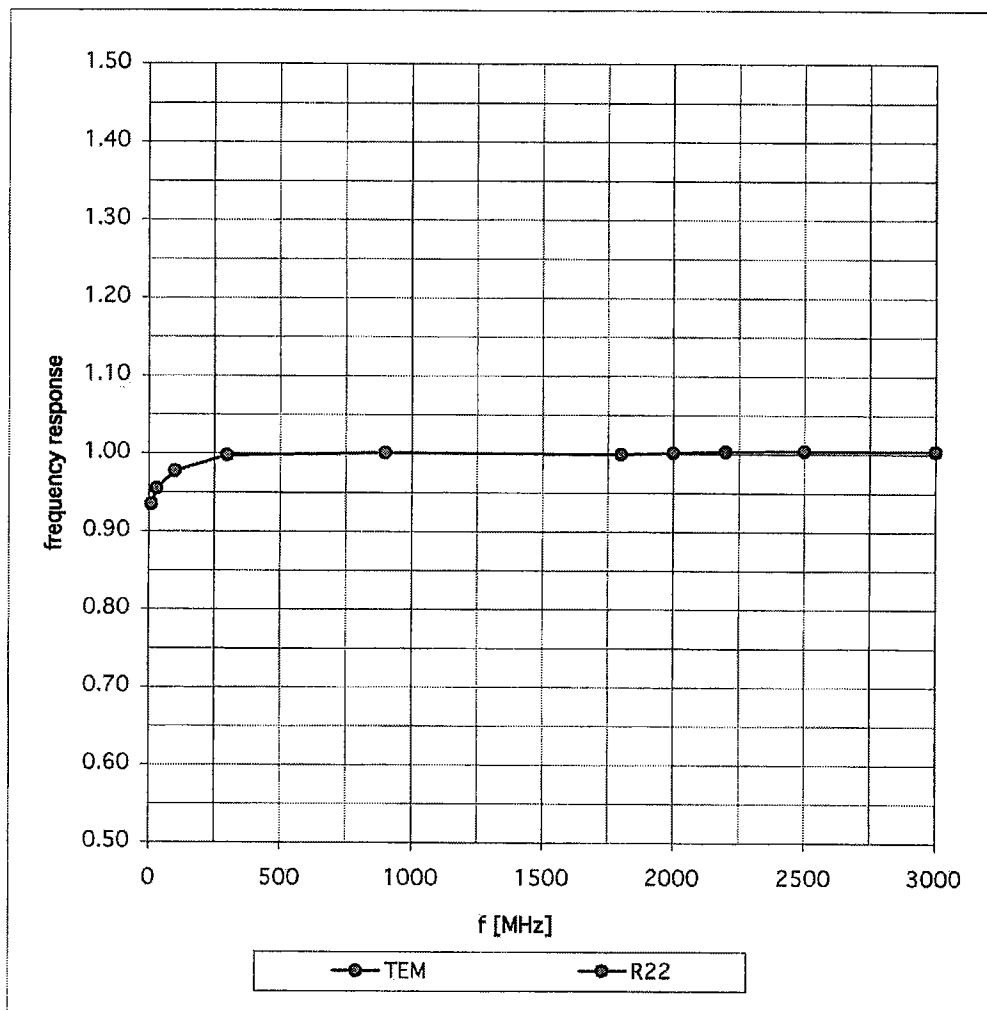


Dynamic Range f(SARhead)

(Waveguide R22)

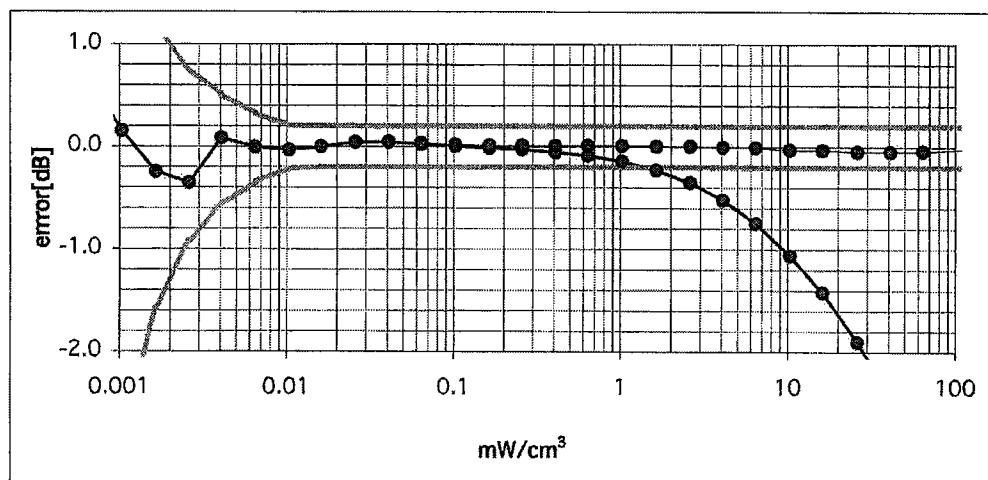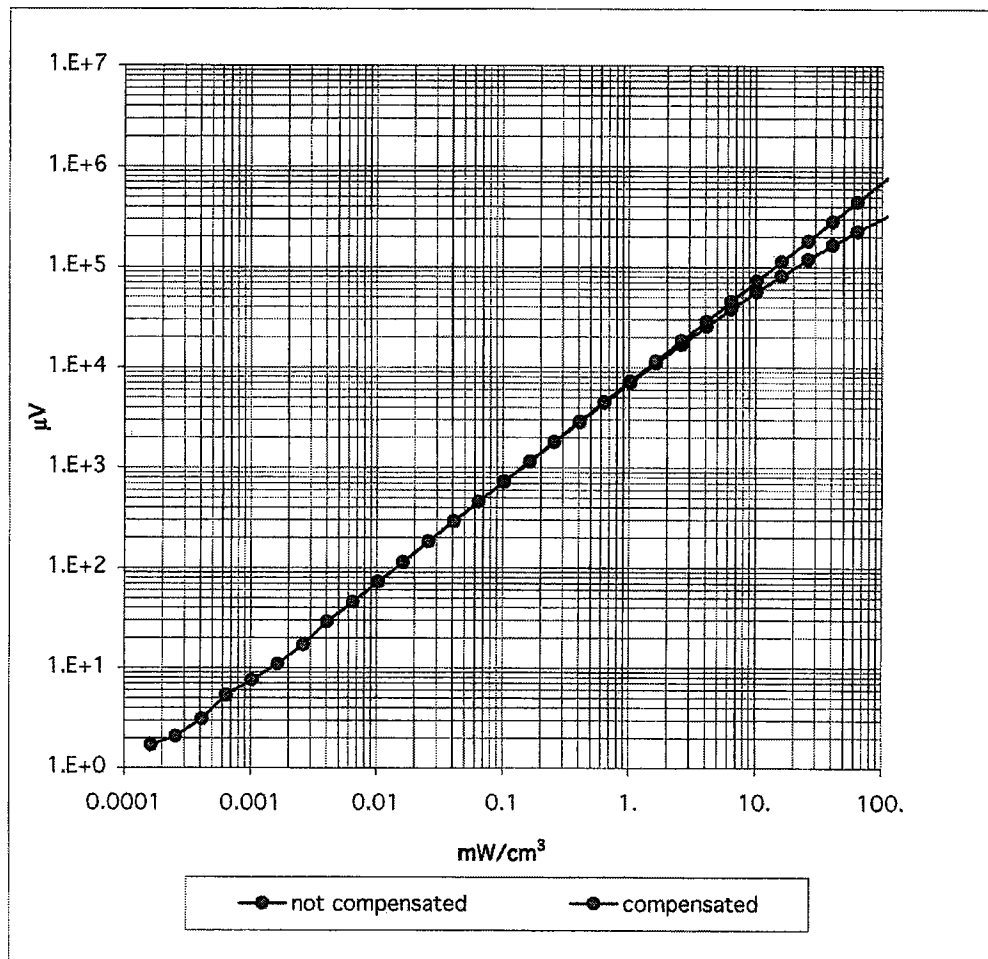

Conversion Factor Assessment

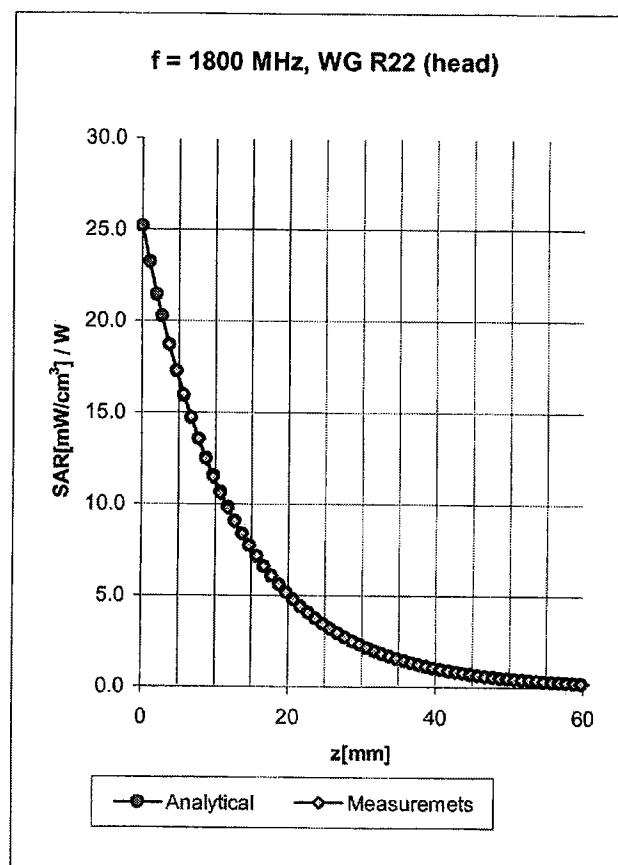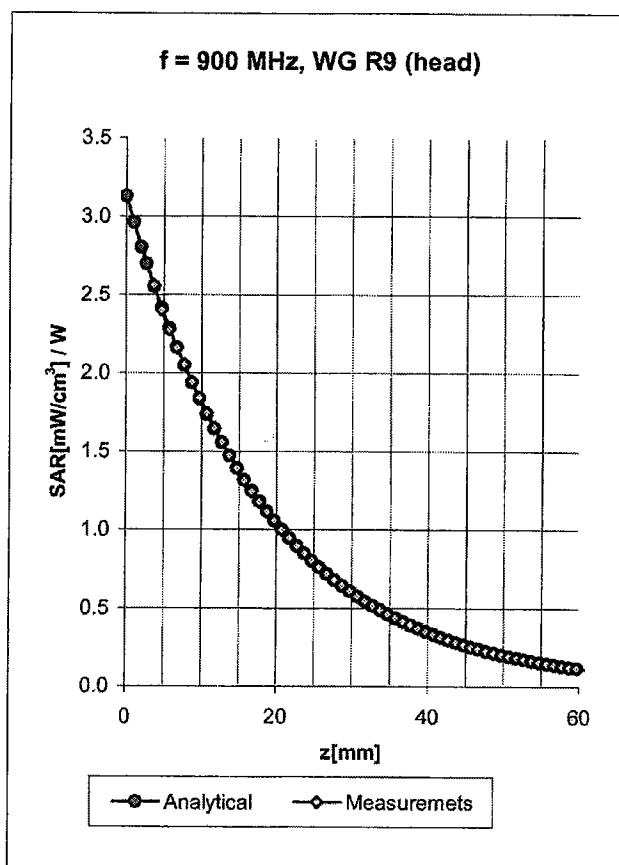

Head 900 MHz $\epsilon_r = 41.5 \pm 5\%$ $\sigma = 0.97 \pm 5\%$ mho/m

Valid for f=800-1000 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	6.3 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	6.3 $\pm 9.5\%$ (k=2)	Alpha	0.40
ConvF Z	6.3 $\pm 9.5\%$ (k=2)	Depth	2.60

Head 1800 MHz $\epsilon_r = 40.0 \pm 5\%$ $\sigma = 1.40 \pm 5\%$ mho/m

Valid for f=1710-1910 MHz with Head Tissue Simulating Liquid according to EN 50361, P1528-200X

ConvF X	5.3 $\pm 9.5\%$ (k=2)	Boundary effect:	
ConvF Y	5.3 $\pm 9.5\%$ (k=2)	Alpha	0.52
ConvF Z	5.3 $\pm 9.5\%$ (k=2)	Depth	2.54

Deviation from Isotropy in HSL

Error (θ, ϕ), $f = 900$ MHz

